


7015A Stainless Steel Deluxe 4-hook I.V. Stand

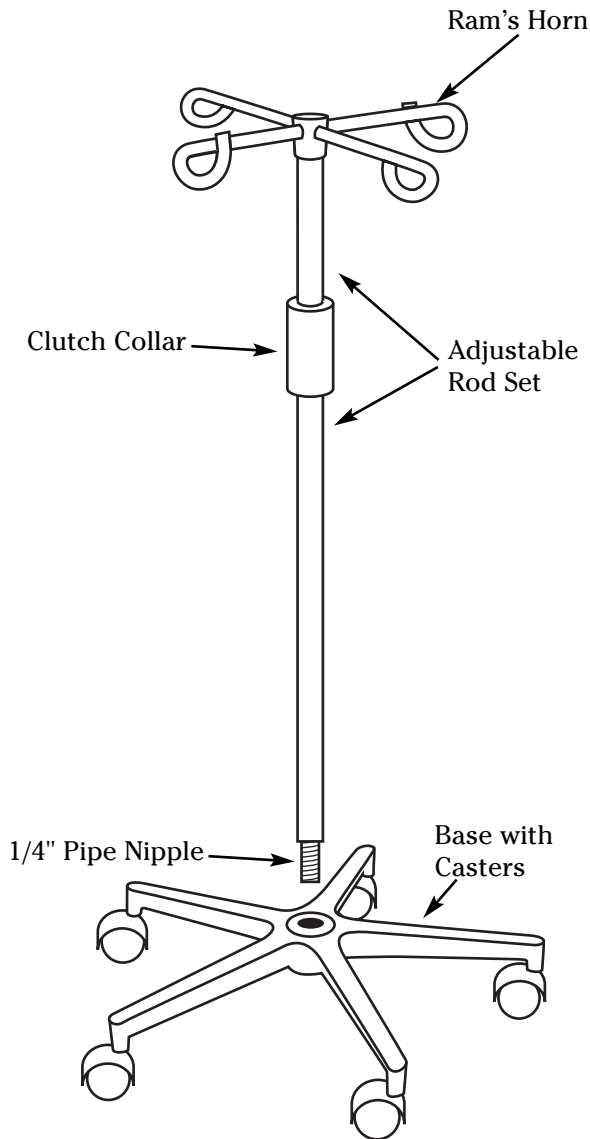

Parts Contents:

Adjustable Rod Set
Clutch Collar
1/4" Pipe Nipple
Base with Casters
1/2" Flat Washer
1/2" External Tooth Lock Washer
Nut
(2) Ram's Horns

Tools Required:

6" Adjustable Wrench

Assembly Steps:

1. Place rod set with 1/4" nipple through base.
2. Place 1/2" flat washer over pipe nipple.
3. Place 1/2" external tooth lock washer over pipe nipple.
4. Secure tightly with nut.
5. Place additional ram's horn in each open hole.

To adjust height, turn collar counter clockwise enough to loosen upper pole. Turn collar clockwise when desired height is obtained.

