


RHYTHM MULTI

Rhythm Multi Mattress System is a medium to high risk system for patients who are at risk of developing Pressure Ulcers according to your sound clinical judgment. Rhythm Multi can also be used for patients who have an existing stage1,2,3 and 4 Pressure Ulcer, in conjunction with your hospitals policy on pressure area management.



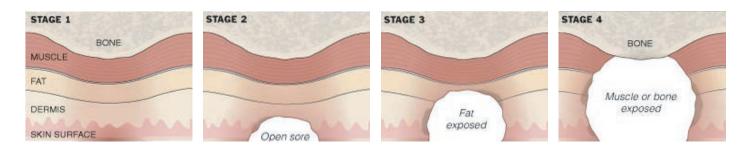

Optional Raised Bolster Cover available for added patient safety. Bolster cover has center access providing easy ingress and egress.

MASTER CONTRO	LUNIT
Dimension:	(L)17.7"x (W)6.8"x (H)10.8"
Weight:	13.2 lbs
Rating Voltage:	110-120 volts
Max Current:	5A
Rating Frequency:	60Hz
Cycle Time:	3~95 mins
	3 in 1 alternation
Min/Max Pressure:	10~61 mmHg +/- 3~4 mmHg
Protection Type:	Class I, Type BF
	Power outage visual and audio alarm
	Static Function / Pulsation
	Auto-firm for fast mattress inflation
	Keypad lock-out function
Air flow rate:	> 1200 liters per minute

MATTRESS REPLACEMENT	
Model No:	Rhythm Multi
Dimensions:	(L)80"x (W)36"x (H)10"
Bariatric Size:	42, 48 inch width
Weight:	28.6 lbs
Weight Capacity:	600, 800, 1000 lbs
Cell Material:	Black Nylon with PU coating
Cell Number:	18 cells
Cover Material:	Blue Nylon with PU coating finish
Base Material:	Black Polyester with PVC coating
	CPR Quick Release System
	Cable Management System on bottom cover
HCPSC Code	E0277

PRIUS HEALTHCARE USA

Creating innovation and difference

T (813) 854-5464 E info@priushc.com

W priushcusa.com

4027 Tampa Road, #3200, OLDSMAR, FL34677

RHYTHM MULTI

Provides you, the caregiver, with peace of mind, and your patient with world class care and comfort



PRIUS HEALTHCARE USA Creating innovation and difference





AT Prius

priushcusa.com

Perfect for



HAPP extra comfort settings for the heels



Low slip nylon under mattress base cover to reduce shearing and friction between mattress and bed frame

Nylon woven fabric with PU coating finish

PRESSURE MAPPING

Rhythm Multi provides constant low pressure by redistributing the patient's pressure over a greater area than a standard mattress.

RHYTHM MULTI

- 1000 lbs
- Nylon PU cell design provides easy cleaning and minimizing the proliferation of bacteria
- 60 seconds of pulsation simulates a message to assist in maximizing a patient's comfort
- Modular structure design to ensures ease of use and cleaning

 - temperature control

2" Foam Safety Mat to prevent bottoming out in case of power loss

-

Exhaust port removes excess heat from pump lowering mattress surface temperature

10 adjustable comfort settings provides different levels of comfort depending on the patient's therapy required

3 different modes of clinical therapy static/ alternation/pulsation



• The Rhythm Multi replacement mattress can accommodate patient weights up to

• True low air loss provides extra comfort for moisture and



provides 3 selectable softness settings in the last 4 air cells to provide protection for patient's fragile heels



Alternation cycle times from 3-95 minutes

Control Panel Lock Out

Auto Firm/Patient Care control



Rhythm Multi is your partner in protection and treatment of pressure ulcer



Cable management keeps electrical cables off the floor



CPR quick release system deflates the system within 20 seconds

