


DynaFlow

DynaFlow Mattress System is a medium to high risk system for patients who are at risk of developing Pressure Ulcers according to your sound clinical judgment. DynaFlow can also be used for patients who have an existing stage 1,2,3 and 4 Pressure Ulcer, in conjunction with your hospitals policy on pressure area management.





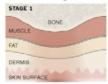









Various top covers are designed and available to satisfy your individual requirement in terms of therapy, prevention and patient comfort.

Mattress Control Unit

- Dimension: (L)12.3"x (W)6.7"x (H)4.9"
- Weight: 14.9 lbs
- Rated Voltage: 110-120 volts
- Max Current: 5A
- Rated Frequency: 60Hz
- Cycle Time: 3~95 mins
- 3 in 1 alternation

- Min/Max Pressure: 15~52 mmHg +/- 4 mmHg
- Static / Pulsation Function

Mattress Replacement

- Model No: DynaFlow
- Dimension: (L)80"x (W)36"x (H)10"
- · Weight: 39.6 lbs
- Weight Capacity: 600lbs
- Cell Material: Blue Nylon with PU Coating

- Protection Type: Class I, Type BF
- 10 Digital Adjustable Comfort Control
- Bottom Out Alarm: Visual / Audio
- Leaving Bed Alarm: Visual / Audio
- Power Failure Alarm: Visual / Audio
- User-Friendly LCD Control Panel
- Fowler Function
- Auto Function

- · Cell Number: 18 cells
- Cover Material: Blue Nylon with PU coating finish.
- Base Material: Black Polyester with PVC Coating
- CPR Quick Release System





CREATING INNOVATION AND DIFFERENCE







provides you, the carer, with peace of mind and your patient, with world class care & comfort

Perfect for patients and caregivers



Low slip base cover is a plus to keep the mattress on the bed

The modular construction of the cells makes for simplistic repair and cleaning



Incorporates with automatic sensing technology for bottoming out adjustment

PRESSURE MAPPING

nedium

DynaFlow provides constant low pressure by redistributing the patient's pressure over a greater area than

a standard mattress.



. Equip with 2" convoluted foam base provides additional support during power failure situation

> PU cell design provides easy cleaning and minimizing the proliferation of bacteria

> > Mattress replacement can accommodate patient weights up to 600 lbs

> > > · Added safety feature for continous monitoring patient's status on bottoming out and leaving bed

> > > > Low Air Loss, Pulsation, 3-1 Alternation. Static and Rotation Therapies in one system

Cable management keeps electrical cables off the floor

DynaFlow

is your partner

in protection

and treatment

of pressure ulcer



release system deflates the system within 20 seconds



Hanging bracket

DynaFlow can detect the pressure of patients' buttocks automatically during "Fowler Mode" and adjust the pressure based on these results.

male 5.5 ft/154.3lbs

Automatically sets optimal pressure based on patient's weight, height and position.

30 seconds of pulsation simulates a massage to assist in maximizing a patients comfort.



The mattress can be adjusted for "Fowler Mode" while acting in conjunction with the bed frame, assisting the medical staff while performing medical procedures.

Multiple lateral rotations provide alternative treatments to patients and medical staffs. Rotation therapy can be set from 3 to 95 minutes.